MAGERI benchmark using reference standard DNA library

Mikhail Shugay
December 1, 2016

```
Load metadata
```

```
df.vmeta <- read.table("hd734_variant_metadata.txt", sep="\t", header=T) # variants observed in HD734 a
df.smeta <- read.table("sample_metadata.txt", sep="\t", header=T) # metadata for amplicon sequencing sa
VCF parsing function
library(stringr)
library(pROC)
## Type 'citation("pROC")' for a citation.
## Attaching package: 'pROC'
## The following objects are masked from 'package:stats':
##
##
       cov, smooth, var
library(ggplot2)
library(ggbeeswarm)
library(dplyr)
## Attaching package: 'dplyr'
## The following objects are masked from 'package:stats':
##
##
       filter, lag
## The following objects are masked from 'package:base':
##
##
       intersect, setdiff, setequal, union
read_vcf <- function(file_name) {</pre>
  .vcf <- read.table(file_name, header = F, sep = "\t", stringsAsFactors = F)</pre>
  colnames(.vcf) <- c("chromosome", "position", "skip1", "from", "to", "qual", "skip2", "info", "skip3"</pre>
  .vcf$skip1 <- NULL</pre>
  .vcf$skip2 <- NULL</pre>
  .vcf$skip3 <- NULL
  .vcf$skip4 <- NULL
```

```
.vcf$qual <- as.integer(.vcf$qual)</pre>
  .vcf <- subset(.vcf, nchar(from) == 1 &</pre>
                  nchar(to) == 1 & !is.na(qual)) # no indels
  .infosplit <- str_split_fixed(.vcf$info, regex("[=;]"), 15)[,c(2, 4, 12, 14)]</pre>
  .vcf$coverage <- as.numeric(.infosplit[,1])</pre>
  .vcf$frequency <- as.numeric(.infosplit[,2])</pre>
  .vcf$bb.a <- as.numeric(.infosplit[,3])</pre>
  .vcf$bb.b <- as.numeric(.infosplit[,4])</pre>
  .vcf$info <- NULL</pre>
  .vcf$count <- as.integer(round(.vcf$coverage * .vcf$frequency))</pre>
  .vcf
}
head(read_vcf("p126.h4_2_ballast_m1.vcf"))
## Warning in read_vcf("p126.h4_2_ballast_m1.vcf"): NAs introduced by coercion
    chromosome position from to qual coverage
                                                frequency
## 1
        chr2 212295704 C A 8
                                        2265 0.0004415011 0.9957058
         chr2 212295705 C A 8
## 2
                                        2265 0.0004415011 0.9957058
## 3
         chr2 212295713 A G 7
                                       2265 0.0004415011 1.0295689
## 4
         ## 5
          chr2 212295732 C T 24
## 6
                                       2264 0.0013250883 2.1314828
##
        bb.b count
## 1 21862.23
## 2 21862.23
## 3 16394.12
                1
## 4 24344.82
                1
## 5 21862.23
                1
## 6 24344.82
                 3
Read samples with HD734 standard DNA and control human DNA, append metadata
library(TailRank) # For betabinom.
## Loading required package: oompaBase
# installing:
# source("https://bioconductor.org/biocLite.R")
# biocLite("Biobase")
# install.packages("TailRank", repos="http://R-Forge.R-project.org")
df <- data.frame()</pre>
read_vcf_with_metadata <- function(file_name, primer_set, replica, ratio, type) {</pre>
```

```
.vcf <- read_vcf(file_name)</pre>
  .vcf <- merge(.vcf, df.vmeta, all.x = type != "standard", all.y = F)</pre>
  .vcf$known.frequency <- .vcf$known.frequency * ratio</pre>
  .vcf\$known.frequency[is.na(.vcf\$known.frequency)] <- 0
  .vcf$primer_set <- primer_set</pre>
  .vcf$replica <- primer_set</pre>
  .vcf$type <- type</pre>
  .vcf <- subset(.vcf, frequency < 0.4 & count > 0) # remove alleles in control
  .vcf
}
for (i in 1:nrow(df.smeta)) {
  df <- with(df.smeta, rbind(df,</pre>
                                read_vcf_with_metadata(paste(prefix[i], "vcf", sep="."),
                                                         primer_set[i],
                                                         replica[i],
                                                         ratio[i],
                                                         type[i])))
}
```

```
## Warning in read_vcf(file_name): NAs introduced by coercion
```

Warning in read_vcf(file_name): NAs introduced by coercion

 $\#df\$qual \leftarrow ifelse(df\$qual \leftarrow 1e-100, 999, -10*log10(df\$qual))$

Group observed variants into tiers

```
df$tier <- cut(df$known.frequency, c(-1, 0, 0.009, 0.02, 1))
levels(df$tier) <- c("error", "0.1%", "1%", "5%")
summary(df$tier)</pre>
```

```
## error 0.1% 1% 5%
## 2468 43 46 12
```

Number of detected variants and sampling


```
df.v.summary <- df.vmeta
df.v.summary$known.frequency <- df.v.summary$known.frequency * 0.1 # 1:9 diluted samples
df.v.summary$ratio <- "1:10"
df.v.summary.2 <- df.vmeta
df.v.summary.2$ratio <- "1:1"
df.v.summary <- rbind(df.v.summary, df.v.summary.2)
df.v.summary$tier <- cut(df.v.summary$known.frequency, c(-1, 0, 0.009, 0.02, 1))
levels(df.v.summary$tier) <- c("error", "0.1%", "1%", "5%")
print(df.v.summary)</pre>
```

```
##
     chromosome position
                                   from to
## 1
          chr7 140453136
                                      A T
                                      A C
## 2
          chr7 140453136
                                      C T
## 3
          chr7 140453137
## 4
          chr7 55241707
                                      G A
## 5
          chr7 55259515
                                      T G
          chr7 55259524
## 6
                                      Т
                                        Α
## 7
          chr7 55249071
                                      C T
## 8
          chr19
                3118942
                                      A T
## 9
          chr2 209113112
                                      СТ
                                      G A
## 10
          chr2 209113113
                                      C T
## 11
          chr12 25378562
## 12
        chr12 25398284
                                      C G
         chr12 25398285
                                      C A
## 13
```

```
С
## 14
           chr12 25398284
## 15
           chr12 25398285
                                            C
                                               G
           chr12 25398285
## 16
                                            C
                                            С
## 17
           chr12
                  25398284
                                               Α
## 18
           chr12
                  25398281
                                            С
## 19
           chr12 25380275
                                           Т
                                               G
## 20
           chr12 25380276
                                           Τ
## 21
            chr1 115256528
                                               Α
## 22
            chr1 115256530
                                           G
                                               Τ
## 23
                                            Т
            chr1 115256529
                                               Α
## 24
            chr1 115256529
                                            G
## 25
            chr3 178936082
                                               Α
## 26
            chr3 178936091
                                            G
                                               Α
## 27
            chr3 178952085
                                               G
## 28
            chr7 55242465 GGGAATTAAGAGAAGC
                                               G
## 29
            chr7 140453136
                                               Т
## 30
                                               С
            chr7 140453136
                                            Α
                                               Т
## 31
            chr7 140453137
                                            С
## 32
                                            G
            chr7 55241707
                                               Α
                                            Τ
## 33
            chr7 55259515
                                               G
## 34
            chr7 55259524
                                            Τ
                                               Α
## 35
            chr7
                  55249071
                                            С
## 36
                                           Α
                                               Τ
           chr19
                   3118942
## 37
            chr2 209113112
                                           C
                                               Τ
## 38
            chr2 209113113
                                           G
                                               Α
## 39
           chr12 25378562
                                           C
                                               Т
## 40
           chr12
                  25398284
                                           С
                                               G
                  25398285
                                            С
## 41
           chr12
                                               Α
                                            C
                                               Τ
## 42
                  25398284
           chr12
                                           C
## 43
           chr12
                  25398285
                                               G
                                           C
## 44
           chr12
                  25398285
                                               Τ
## 45
           chr12
                  25398284
                                           C
                                               Α
                                            С
                                               Т
## 46
           chr12 25398281
## 47
                  25380275
                                           Τ
                                               G
           chr12
                                            Τ
## 48
           chr12 25380276
## 49
            chr1 115256528
                                           Т
                                               Α
## 50
            chr1 115256530
                                            G
## 51
            chr1 115256529
                                           Τ
                                               Α
## 52
            chr1 115256529
                                            Т
                                               C
## 53
                                            G
            chr3 178936082
## 54
            chr3 178936091
                                               Α
## 55
            chr3 178952085
                                            Α
                                               G
## 56
            chr7 55242465 GGGAATTAAGAGAAGC
                                              G
##
                                     id known.frequency
                                                          reference ratio tier
                 BRAF_V600E, BRAF_V600K
## 1
                                                  0.0080
                                                           BRAF_E15 1:10 0.1%
                 BRAF_V600G, BRAF_V600R
                                                           BRAF_E15
## 2
                                                  0.0010
                                                                     1:10 0.1%
      BRAF_V600M, BRAF_V600K, BRAF_V600R
## 3
                                                  0.0010
                                                           BRAF_E15
                                                                     1:10 0.1%
## 4
                                                           EGFR_E18
                             EGFR_G719S
                                                  0.0167
                                                                     1:10
## 5
                             EGFR_L858R
                                                  0.0010
                                                           EGFR_E21
                                                                     1:10 0.1%
## 6
                             EGFR_L861Q
                                                  0.0010
                                                           EGFR_E21
                                                                      1:10 0.1%
                             EGFR_T790M
## 7
                                                           EGFR_E20
                                                  0.0010
                                                                     1:10 0.1%
## 8
                            GNA11_Q209L
                                                  0.0013
                                                           GNA11_E5 1:10 0.1%
## 9
                             IDH1_R132C
                                                  0.0013
                                                            IDH1_E4 1:10 0.1%
## 10
                             IDH1_R132H
                                                  0.0013
                                                            IDH1 E4 1:10 0.1%
```

```
## 11
                              KRAS A146T
                                                    0.0013
                                                              KRAS E4
                                                                       1:10 0.1%
## 12
                               KRAS_G12A
                                                    0.0013
                                                              KRAS E2
                                                                        1:10 0.1%
## 13
                               KRAS G12C
                                                    0.0013
                                                              KRAS E2
                                                                        1:10 0.1%
                               KRAS_G12D
                                                              KRAS_E2
## 14
                                                    0.0013
                                                                        1:10 0.1%
## 15
                               KRAS_G12R
                                                    0.0013
                                                              KRAS_E2
                                                                        1:10 0.1%
## 16
                               KRAS G12S
                                                    0.0013
                                                              KRAS E2
                                                                       1:10 0.1%
## 17
                               KRAS G12V
                                                    0.0013
                                                              KRAS E2
                                                                        1:10 0.1%
## 18
                               KRAS_G13D
                                                    0.0250
                                                              KRAS_E2
                                                                        1:10
## 19
                               KRAS_Q61H
                                                    0.0013
                                                              KRAS_E3
                                                                        1:10 0.1%
## 20
                               KRAS_Q61L
                                                    0.0013
                                                              KRAS_E3
                                                                        1:10 0.1%
## 21
                               NRAS_Q61H
                                                    0.0013
                                                              NRAS_E3
                                                                        1:10 0.1%
## 22
                               NRAS_Q61K
                                                    0.0013
                                                              NRAS_E3
                                                                        1:10 0.1%
## 23
                               NRAS_Q61L
                                                    0.0013
                                                              NRAS_E3
                                                                        1:10 0.1%
                                                                        1:10 0.1%
## 24
                               NRAS_Q61R
                                                    0.0013
                                                              NRAS_E3
## 25
                            PIK3CA_E542K
                                                    0.0013
                                                            PIK3CA_E9
                                                                        1:10 0.1%
## 26
                           PIK3CA_E545K
                                                    0.0013
                                                            PIK3CA_E9
                                                                        1:10 0.1%
## 27
                                                    0.0300 PIK3CA_E20
                                                                                5%
                           PIK3CA_H1047R
                                                                        1:10
## 28
                        EGFR dE746-A750
                                                    0.0010
                                                             EGFR E19
                                                                        1:10 0.1%
## 29
                  BRAF_V600E, BRAF_V600K
                                                    0.0800
                                                             BRAF_E15
                                                                         1:1
                                                                                5%
## 30
                  BRAF V600G, BRAF V600R
                                                    0.0100
                                                             BRAF E15
                                                                         1:1
                                                                                1%
## 31
      BRAF_V600M, BRAF_V600K, BRAF_V600R
                                                    0.0100
                                                             BRAF_E15
                                                                         1:1
                                                                                1%
                              EGFR G719S
                                                    0.1670
                                                             EGFR E18
                                                                         1:1
                                                                                5%
## 33
                              EGFR_L858R
                                                    0.0100
                                                             EGFR_E21
                                                                         1:1
                                                                                1%
## 34
                              EGFR L861Q
                                                    0.0100
                                                             EGFR E21
                                                                         1:1
                                                                                1%
## 35
                              EGFR T790M
                                                    0.0100
                                                             EGFR E20
                                                                         1:1
## 36
                             GNA11_Q209L
                                                    0.0130
                                                             {\tt GNA11\_E5}
                                                                         1:1
                                                                                1%
## 37
                              IDH1_R132C
                                                    0.0130
                                                              IDH1_E4
                                                                         1:1
                                                                                1%
## 38
                              IDH1_R132H
                                                    0.0130
                                                              IDH1_E4
                                                                         1:1
                                                                                1%
## 39
                              KRAS_A146T
                                                    0.0130
                                                              KRAS_E4
                                                                         1:1
                                                                                1%
## 40
                               KRAS_G12A
                                                              KRAS_E2
                                                                         1:1
                                                    0.0130
                                                                                1%
## 41
                               KRAS_G12C
                                                    0.0130
                                                              KRAS_E2
                                                                         1:1
                                                                                1%
## 42
                               KRAS_G12D
                                                    0.0130
                                                              KRAS_E2
                                                                         1:1
                                                                                1%
## 43
                               KRAS_G12R
                                                    0.0130
                                                              KRAS_E2
                                                                         1:1
                                                                                1%
## 44
                               KRAS_G12S
                                                    0.0130
                                                              KRAS_E2
                                                                         1:1
                                                                                1%
## 45
                                                    0.0130
                                                              KRAS E2
                                                                         1:1
                                                                                1%
                               KRAS_G12V
## 46
                               KRAS_G13D
                                                    0.2500
                                                              KRAS_E2
                                                                         1:1
                                                                                5%
## 47
                               KRAS Q61H
                                                    0.0130
                                                              KRAS E3
                                                                                1%
## 48
                               KRAS_Q61L
                                                              KRAS_E3
                                                                         1:1
                                                                                1%
                                                    0.0130
                                                    0.0130
                                                                                1%
## 49
                               NRAS_Q61H
                                                              NRAS E3
                                                                         1:1
## 50
                               NRAS_Q61K
                                                    0.0130
                                                              NRAS_E3
                                                                         1:1
## 51
                               NRAS Q61L
                                                    0.0130
                                                              NRAS E3
                                                                         1:1
                                                                                1%
## 52
                               NRAS_Q61R
                                                    0.0130
                                                              NRAS E3
                                                                         1:1
                                                                                1%
## 53
                           PIK3CA_E542K
                                                    0.0130
                                                            PIK3CA E9
                                                                         1:1
                                                                                1%
## 54
                                                                                1%
                            PIK3CA_E545K
                                                    0.0130
                                                            PIK3CA_E9
                                                                         1:1
## 55
                           PIK3CA_H1047R
                                                    0.3000 PIK3CA_E20
                                                                         1:1
                                                                                5%
## 56
                                                   0.0100
                                                                         1:1
                                                                                1%
                         EGFR_dE746-A750
                                                             EGFR_E19
df.v.summary <- df.v.summary %>%
  group_by(tier) %>%
  summarise(total = n())
df.v.summary$total <- 2 * df.v.summary$total # we have 2 independent experiments
df.v.summary <- merge(df.v.summary,</pre>
                       df %>% group_by(tier) %>%
```

```
summarise(found = n(),
                                  detect.prob.m = mean(1-exp(-frequency * coverage)),
                                   detect.prob.sd = sd(1-exp(-frequency * coverage))))
df.v.summary$detection.prob.lb <- with(df.v.summary,</pre>
                                        detect.prob.m - 1.96*detect.prob.sd - 5e-3)
df.v.summary$detection.prob.ub <- with(df.v.summary,</pre>
                                        detect.prob.m + 1.96*detect.prob.sd)
df.v.summary$detection.prob.ub <- ifelse(df.v.summary$detection.prob.ub > 1,
                                          df.v.summary$detection.prob.ub )
library(scales)
fig3 <- ggplot(df.v.summary, aes(x=tier)) +</pre>
  geom_bar(stat="identity", aes(y=detection.prob.ub), fill = "grey70", width=0.5)+
  geom_bar(stat="identity", aes(y=detection.prob.lb), fill = "white")+
  geom_point(aes(y = found/total), size = 5, color="red", shape="+") +
  geom_text(aes(y = found/total - 0.05, label = paste(found, total, sep="/"))) +
  scale_y_continuous("Detection rate", labels=percent) +
 xlab("Variant tier") +
 theme bw() +
 theme(panel.grid.major.x = element_blank(),
        panel.grid.major.y = element_blank())
save(file = "../figures/fig3.Rda", fig3)
fig3
```

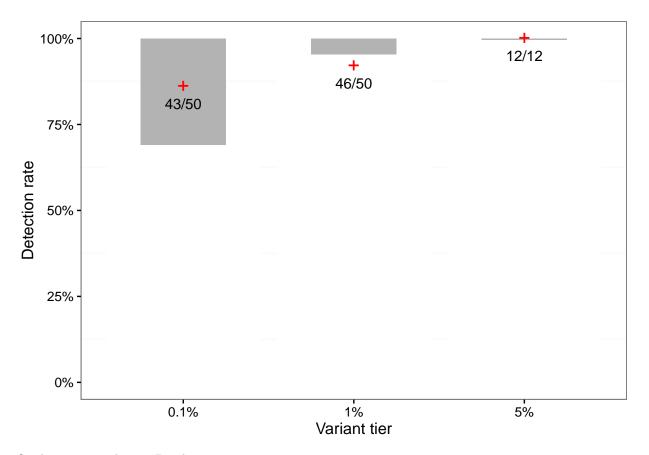


Quality score and error P-values

```
df.e <- df %>%
  filter(type == "blank") %>%
  mutate(true.p.value = -10 * log10(1 - (rank(frequency) - 0.5) / n()))

fig4 <- ggplot(df.e, aes(true.p.value, qual)) +
  geom_abline(slope = 1, intercept = 0, linetype = "dashed") +
  geom_smooth(color="red", fill="grey", alpha=0.5) +
  scale_x_continuous("-10 log10 P-value", limits = c(0, 40)) +
  scale_y_continuous("Q score", limits = c(0, 40)) +
  theme_bw()

save(file = "../figures/fig4.Rda", fig4)</pre>
```


Warning: Removed 1 rows containing non-finite values (stat_smooth).

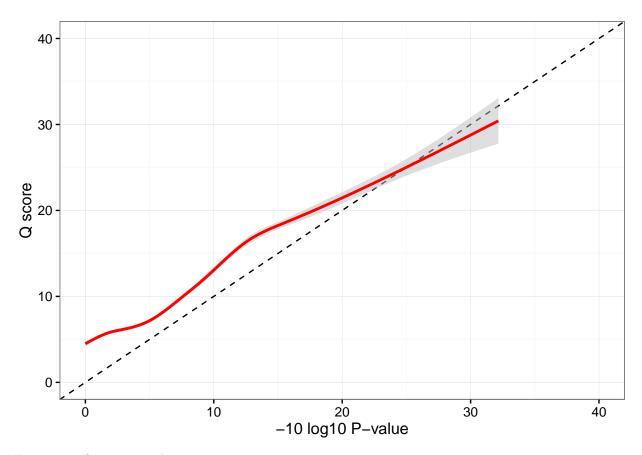

Frequency of variants and errors

```
fig5 <- ggplot(df, aes(x=frequency, color = tier)) +
  geom_density(aes(fill=tier), alpha=0.7) +
  ylab("Density") +
  scale_x_log10("Frequency") +
  scale_color_brewer("Variant tier", palette = "Set1") +
  scale_fill_brewer("Variant tier", palette = "Set1") +
  theme_bw()


save(file = "../figures/fig5.Rda", fig5)</pre>
```

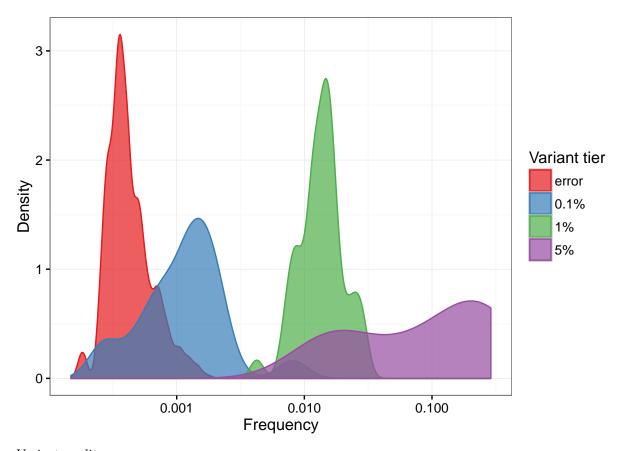

Variant quality scores

```
fig6 <- ggplot(df,aes(tier, qual, fill=tier)) +
   geom_hline(yintercept = 13.0103, linetype ="dotted") +
   geom_hline(yintercept = 20.0103, linetype ="dashed") +
   geom_quasirandom(varwidth = T, shape=21, color="grey10") +
   annotate("text", label = "p < 0.05", x=4.3,y=13.0103+2) +
   annotate("text", label = "p < 0.01", x=4.3,y=20.0103+4) +
   scale_fill_brewer("", palette = "Set1", guide=F) +
   xlab("Variant tier") + ylab("Q score") +
   theme_bw()

save(file = "../figures/fig6.Rda", fig6)</pre>
```

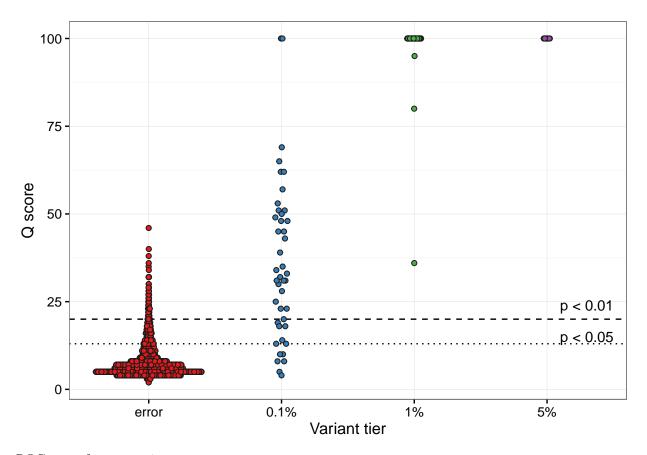

ROC curve for rare variants

Area under the curve: 0.9277 ## 95% CI: 0.8727-0.9827 (DeLong)

```
df.1 <- subset(df, tier %in% c("error", "0.1%"))</pre>
df.1$type <- ifelse(df.1$tier == "error", 0, 1)</pre>
make_roc <- function(rocobj, type) {</pre>
  .df.roc \leftarrow data.frame(spec = seq(0, 1, 0.01))
  sens.ci <- ci.se(rocobj, specificities = .df.roc$spec)</pre>
  .df.roc$sens.lo <- sens.ci[,1]</pre>
  .df.roc$sens.me <- sens.ci[,2]</pre>
  .df.roc$sens.hi <- sens.ci[,3]</pre>
  .df.roc$type <- type</pre>
  .df.roc
}
rocobj <- roc(type ~ qual, df.1, ci=T)</pre>
print(rocobj)
##
## Call:
## roc.formula(formula = type ~ qual, data = df.1, ci = T)
## Data: qual in 2468 controls (type 0) < 43 cases (type 1).
```

```
df.roc <- make_roc(rocobj, "Q score")</pre>
rocobj <- roc(type ~ frequency, df.1, ci=T)</pre>
print(rocobj)
##
## Call:
## roc.formula(formula = type ~ frequency, data = df.1, ci = T)
##
## Data: frequency in 2468 controls (type 0) < 43 cases (type 1).
## Area under the curve: 0.8637
## 95% CI: 0.7839-0.9435 (DeLong)
df.roc <- rbind(df.roc, make_roc(rocobj, "Frequency"))</pre>
fig7 <- ggplot(df.roc, aes(x=spec)) +</pre>
  geom_ribbon(aes(ymin=sens.lo, ymax=sens.hi, group=type), fill="grey", alpha=0.5) +
  geom_abline(slope = 1, intercept = 1, linetype = "dashed") +
  geom_line(aes(y=sens.me, color = type)) +
  scale_x_reverse("Specificity") +
  scale_y_continuous("Sensitivity", limits=c(0,1)) +
  scale_color_brewer("Threshold", palette = "Set1") +
  theme bw()
save(file = "../figures/fig7.Rda", fig7)
fig7
```

